

ENFit®  
Compatible

**RELiZORB™** 

(iMMOBILIZED LIPASE) CARTRIDGE


## How to Connect ENFit-Compatible RELiZORB to a Non-ENFit Extension Set




If your extension set  
looks like this,  
follow these important  
instructions.

# HOW TO CONNECT ENFit-COMPATIBLE RELiZORB TO A NON-ENFit EXTENSION SET

Which of these does your extension set look like?

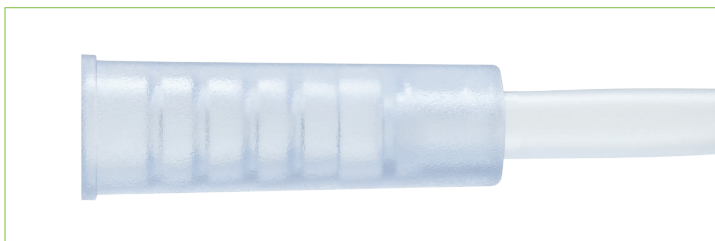

OR

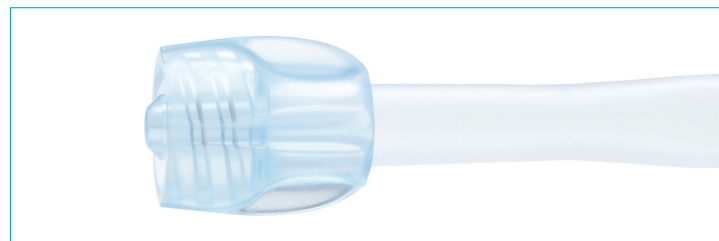

Not ENFit Compatible

ENFit Extension Set

You'll need to remove the white Transition Adapter from the tubing set and move it to the RELiZORB outlet.

Your system is already compatible! Discard the white Transition Adapter and connect the purple fitting from the tubing set to the clear RELiZORB inlet, and then connect RELiZORB to the ENFit extension set.



**STEP 1:** Remove the white Transition Adapter from the tubing set and connect to the outlet of RELiZORB



ENFit completely connected



**STEP 2:** Connect RELiZORB with the transition connector to your extension set.

Questions? Call 1-844-RELiZORB (1-844-735-4967)

Please see Instructions  
for Use for full  
safety information.

**RELiZORB™**   
(IMMOBILIZED LIPASE) CARTRIDGE



## What is ENFit? A New Standard in Enteral Connectors

ENFit is an enteral-specific connector designed to help reduce the risk of enteral tube feeding misconnections and improve patient safety. The ENFit connector helps ensure that connectors do not connect to ports other than ENFit compatible tube feeding connectors.

[www.relizorb.com](http://www.relizorb.com)

1-844-RELiZORB (1-844-735-4967)

